




## MADE IN GERMANY

Wassermann is a medium-sized, owner-operated company with head office in Hamburg. We design and produce high-quality technical devices in Germany. Our main users are dental laboratories, dental clinics and various industrial areas worldwide.

### Why a Wassermann product?

Wassermann products are manufactured from high-quality materials and the proportion of plastics used is intentionally kept as low as possible. Our devices are known for their suitability for daily use, functionality and longevity. As a result, the users place high levels of trust in the Wassermann brand. We feel it is our duty to support our customers on a daily basis when it comes to safety, service and productivity. Technological perfection, optimum quality and especially sustainability are the benchmarks we set. Thanks to the combination of these values, the family-owned company, since its establishment in 1927, has gained worldwide recognition as a manufacturer of high-quality products.

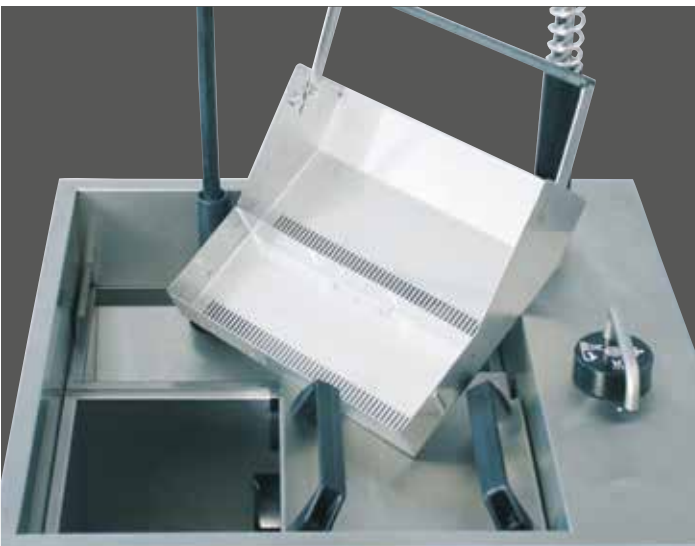

## Wapo-Ex 12 II



## Boil Out Unit



Further information are available on our website at [www.wassermann.hamburg](http://www.wassermann.hamburg)

Your Dental Dealer \_\_\_\_\_



Wassermann Dental-Maschinen GmbH  
Rudorffweg 15-17 · 21031 Hamburg · Germany  
Phone: +49 (0)40 730 926-0 · Fax: +49 (0)40 730 37 24  
info@wassermann-dental.com · www.wassermann.hamburg

# Boil Out Unit **Wapo-Ex 12 II**

The boil out unit **Wapo-Ex 12 II** stands out thanks to its efficient applications planning, programmable timer and automatic wax removal. This free standing unit with integrated spray gallows scalds up to 12 flask halves simultaneously and fully automatically. A touch screen offers diverse functions and reproducible procedures by simple programming. An elaborate boil out unit manufactured using high quality materials and components and designed for durability.

Integrated spray gallows  
for manual scalding of small parts

Selector lever for different functions as  
automatic scalding, manual scalding  
(hand spray) and pump down

Touch screen for  
fast and easy operation

Energy saving  
thermal insulation

High quality construction,  
most parts made of stainless steel



## Major Features

- Free standing unit for scalding up to 12 flask halves
- Efficient usage thanks to its programmable timer and automatic wax removal
- Safety features with visual and audible signals, e.g. electronic water level gauge, leakage signal, inlet limitation
- Energy saving because of the intelligent heating control
- **Included in delivery:**
  - Spray gallows
  - Service set (1 L lime remover and 2 L wax remover)
  - Supply hose with water inlet valve and gasket, 2 m
  - Discharge hose with hose clip Ø 12-20 x 9, 2.5 m
  - Spray arm tube with seal
  - 2 flask baskets
  - 2 covers

### Accessories

Service set (1 L lime remover and 2 L wax remover)

### Item no.

170330

### Wapo-Ex 12 II

Voltage

220-240 V / 50/60 Hz

Power consumption

12 A

Output

max. 3300 W

W x H (without spray gallows) x D

650 x 900 x 620 mm

W x H (with spray gallows) x D

650 x 1600 x 620 mm

Weight

87 kg

Sound level

≤ 70 dB (A)

Water capacity

approx. 50 l

Water feed rate

4.5 l / min

Water drain height

470 mm